

VIE SCOPE™

EMS



**Designed for
EMS and
Combat
Medics**



To schedule a demo call: 425-577-2713

or visit: www.adroitsurgical.com

Enhanced Line of Sight View

Allows for an excellent straight line of sight view with 360 degree maximal illumination

Easy to Use

Provides the ability to intubate the patient when awake in trauma situations

Patented Lifesaving Technology

Gives the user the best chance of intubating patients in any situation, in or outside of a hospital

A Simple and Efficient Device to Secure the Difficult Airway

Insert scope and identify glottis opening between vocal cords


Pass the bougie under direct vision between the vocal cords into the trachea

Remove the scope leaving the bougie in place

Pass the ETT over the bougie into the trachea

Remove the bougie

Ventilate

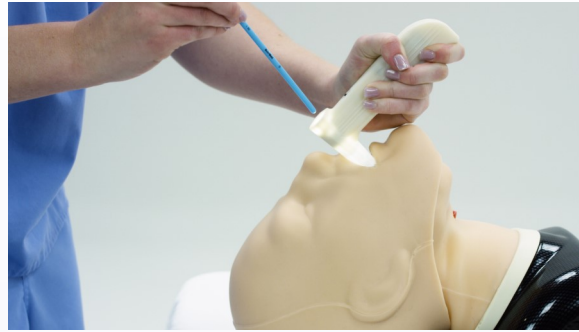

VIE SCOPE™

The Vie Scope is a self-contained, battery powered, disposable scope that takes advantage of a closed circular tube with a beveled end to visualize the vocal cords better than all other methods using the naked eye.

The ring of LED light is transmitted through the side wall of the tube down the entire length of the tube to give the user maximal illumination of the target tissue with minimal chance of light obstruction by secretions or blood.

You could save a life with a Vie Scope.



A Simple and Efficient Device to Secure the Difficult Airway