

VIE SCOPE®

A simple and efficient device to secure the Airway



ENHANCED LINE OF SIGHT VIEW

Allows for an excellent straight line of sight view with 360 degree maximal illumination

EASY TO USE

Provides the ability to intubate the patient when awake in trauma situations

PATENTED LIFESAVING TECHNOLOGY

Gives the users the best chance of intubating patients in any situation, in or outside of a hospital

**To schedule a demo call: 1-877-7ADROIT
or visit: www.adroitsurgical.com**

Adroit Surgical LLC -13904 Quailbrook Drive Oklahoma City - OK 73134

www.adroitsurgical.com



VIE SCOPE®

A simple and efficient device to secure the Airway



Insert scope and identify glottis opening
between vocal cords


Pass the bougie under direct vision
between the vocal cords into the trachea

Remove the scope leave the bougie
in place

Pass the ETT over the bougie into the
trachea

Remove the bougie

Ventilate

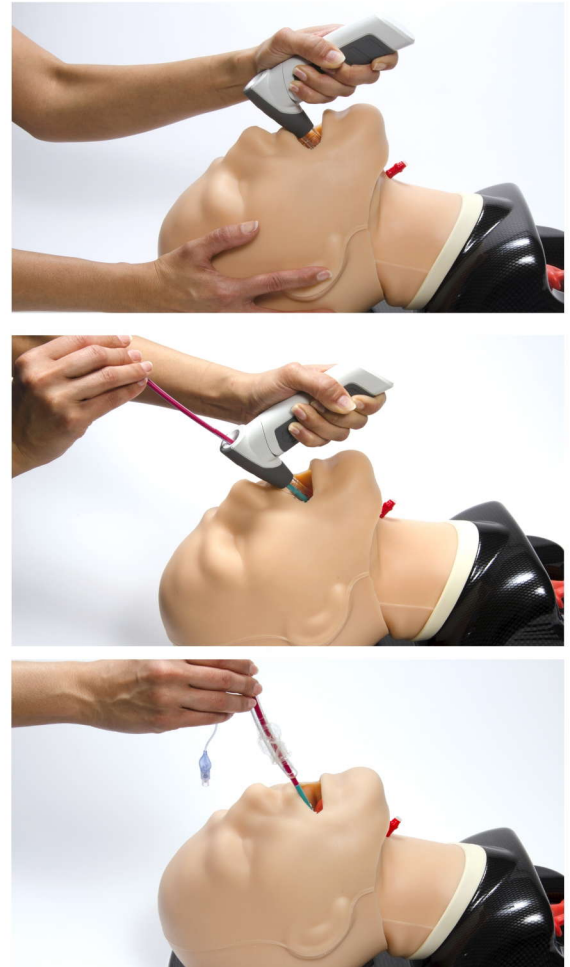

The Vie Scope is a self-contained, battery powered, disposable scope that takes advantage of a closed circular tube with a beveled end to visualize the vocal cords better than all other methods using the naked eye.

The ring of LED lights is transmitted through the side wall of the tube down the entire length of the tube to give the user maximal illumination of the target tissue with minimal chance of light obstruction by secretions or blood.

You could save a life with a Vie Scope.

Adroit Surgical LLC -13904 Quailbrook Drive Oklahoma City - OK 73134

www.adroitsurgical.com

